ADVERTENCIA

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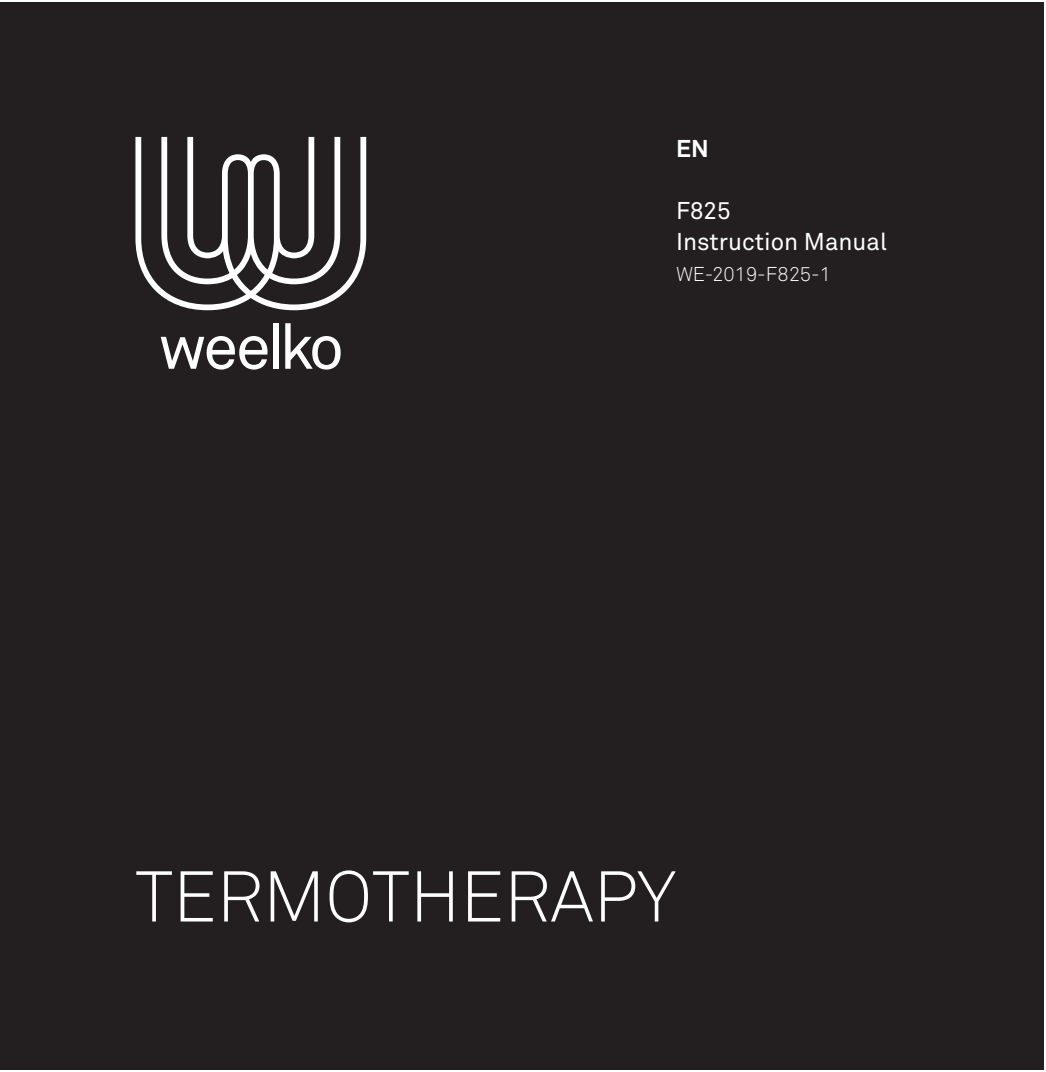
F825
Instruction Manual
WE-2019-F825-1

WARNING

- Never under any circumstances attempt to open or inspect the internal components or accessories of the machine. If an inspection becomes necessary, please contact your supplier or a suitably competent technician.
- Never use the unit in humid or poorly ventilated environments, such as wet rooms, bathrooms, etc. to prevent damage to the electrical components of your unit.
- Never use the unit outdoors; exposure to the elements can cause damage.
- If the power cord was destroyed, please ask a professional to service and replace it to avoid hazards, or contact the dealer for service.

EC-DECLARATION OF COMPLIANCE

This product meets the following requirements:

1. Conditions set out in the Low Voltage Directive: 2014/35/EC
2. Conditions set out in the Electromagnetic Compatibility Directive: 2014/30/EC

RECYCLING

 Whenever any of our electrical products or their presentation packaging bear this symbol, it means that these products should not be treated as conventional household waste in Europe. To ensure proper treatment of this waste, dispose of it in accordance with local laws or as required for the disposal of electrical equipment. This initiative will help preserve save natural resources and improve standards of environmental protection in regard to the treatment of electrical waste.

WELCOME!

Thank you very much for purchase. Please read the instructions carefully before operating the instrument. You can also ask a professional for guidance to make sure you are using the instrument correctly. We will accept no responsibility




for any accidents that occur due to incorrect operation.

The company reserves the right to amend the details of the instruments with no prior notice. If any mistakes are found in these instructions, please feel free to send us an e-mail informing us so that we can correct them.



PACKAGE CONTENTS

- These components are subject to change without prior notice.
- Appearance of product may differ from the image.

Nr.	PRODUCT	REF.	QTY.
1	THERMAL BLANKET	40002	1
2	MAIN UNIT	-	1
3	FEEDER CABLE	-	1

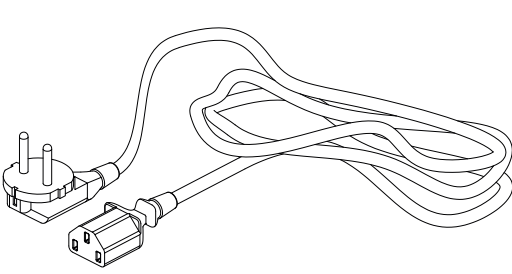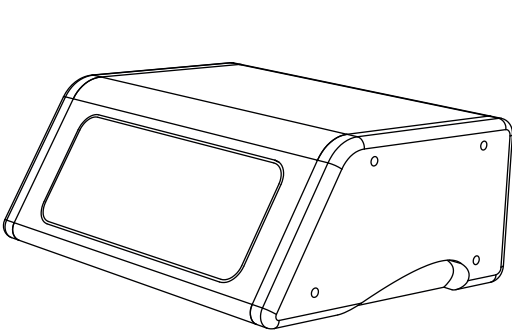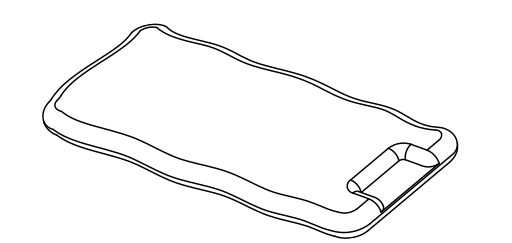FRONT PANEL/REAR PANEL

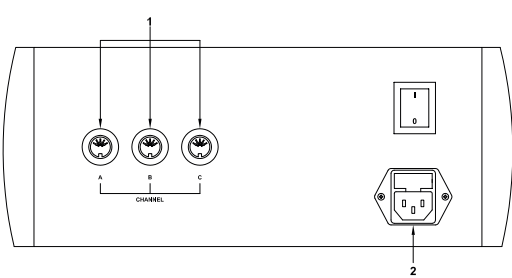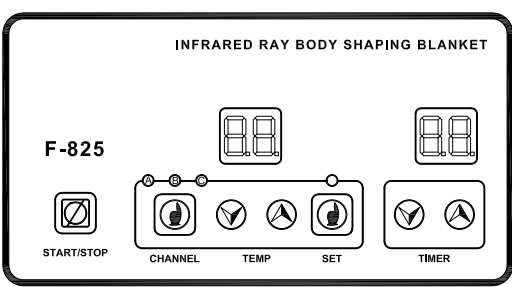1. Channel A/B/C

2. Power line input

INSTRUCTIONS OF USE

1. Connect the power supply via the rear panel, indicated in the sketch as “POWER LINE INPUT”.
2. Connect the accessories properly to the front/rear panel as indicated in the sketch.
3. Turn on the “POWER” button.
4. The unit is in standby. The “TEMP” button is locked, but “TIMER”, “START/STOP”, “SET” and “CHANNEL” are unlocked.

5. You can press the “CHANNEL” button to select one of the 3 different channels. The temperature display shows the actual temperature for each channel.

6. Press the “TIMER” button to set the time from 1 to 60 minutes.

7. Press the “SET” button: the instrument enters setting mode and the “SET” LED comes on. Now you can adjust the operating temperature by pressing “TEMP” in each channel. The default temperature for the three channels is 50°C. The adjustment range is from 20 to 50°C.

8. Press the “SET” button to return to standby. The “SET” LED turns off.

9. Press the “START/STOP” button.

10. The unit is operational. The various channels and current temperatures are automatically displayed. The indicator of each channel is red if this channel is in heating state. The indicator will change to green if the channel is in heat preservation state. The “TEMP” and “CHANNEL” buttons are locked, but “TIME”, “START/STOP” and “SET” are unlocked.

11. While operational, if you press “SET”, you can adjust the operating temperature but the unit will stop. To resume the treatment, you must press “START/STOP”.

12. When preset temperature is reached, the unit enters heat preservation mode.

13. When preset time is reached, the unit emits a buzz sound and stops.

CAUTIONS

- If the surface of the blanket is broken or aging, please replace it with a new blanket to prevent accidental skin burns or electric shocks.
- The heating blanket can generate light magnetic radiation to the human body due to its high power. Please note that pregnant women, children, elderly people and people with sensitive skin should avoid using this instrument. Please drink plenty of water before treatment to avoid dehydration and heatstroke.

FAQ

If the buzzer sounds continuously:

- The temperature has reached the maximum allowed temperature (55°C).
- Turn the power off immediately.
- You can reboot the instrument after 30 minutes. Please make sure nobody stays in the heating blanket.
- If this situation has happened before, please contact your dealer.

If temperature display shows “NO”:

- Check all connections.
- If this situation has happened before, please contact your dealer.

MAINTENANCE

- The leather surface of this blanket is made of PVC material. Please keep it dry and clean after use. If the leather is stained with weight-loss cream or other oily substances, use a neutral cleaning liquid to wipe them off. Never use acid or alkaline organic liquids for cleaning.
- Please cut off the power if the instrument will not be used for a long time. Store the housing and accessories separately in a dry environment.

ACCESSORIES

REF.	DESCRIPTION
40002	THERMAL BLANKET

TECHNICAL SHEET

Weights and sizes	Packing size (cm)	67x30,5x40
	Gross weight (kg)	10,5
	Product size (cm)	24x17/24x13
	Blanket size (cm)	180x80
Working characteristics	Working power	450 W
	Temperature range	20 - 50°C
	Maximum temperature	55°C
	Preset working time	30 min
Electrical characteristics	Rated voltage	220V ~ 240V
	Frequency	50Hz/60Hz
	Rated power	580 W

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